

**SIMPLE...
RELIABLE...
VALUE-PRICED**

STERI-VAC CENTRAL EVACUATION UNIT



Includes spare motor

Steri-Vac is designed for affordability, reliability, and easy installation. Its waterless operation is highly economical.

Ultra-Quiet

Steri-Vac's acoustical housing means quiet operation. This model includes a collection canister to separate solids, permitting recovery of lost crowns and inlays.

Uses No Water

This environmentally friendly system helps to conserve a precious resource. Injection pump systems may waste as many as 600 gallons per day.

Simple, Straightforward Installation

Since running water is not used in this system, installation is simple — no water pump to install. Systems are available for 115V or 220-230V to simplify electrical requirements.

No Down Time

A spare motor is supplied, which can easily be changed *without tools* in minutes — no need to call a dealer service technician.

High-Volume Evacuation

This system effectively serves multiple operatories. Update your present system or use Steri-Vac as a replacement power source for your current central power unit.


Specifications

Heavy-duty Motor —
1 1/2 h.p. - 11 amps - 125 cfm
Standard Motor —
7/8 h.p. - 7.9 amps - 95 cfm
Overall Dimensions —
16" wide x 16" deep x 18" high
Shipping weight — 70 lbs.

STERI-DENT

FROM
ALLIED
DIAGNOSTIC IMAGING RESOURCES, INC.
MAKERS OF
REDI-CHEM

STERI-DENT - Oral evacuation mouthpieces and accessories

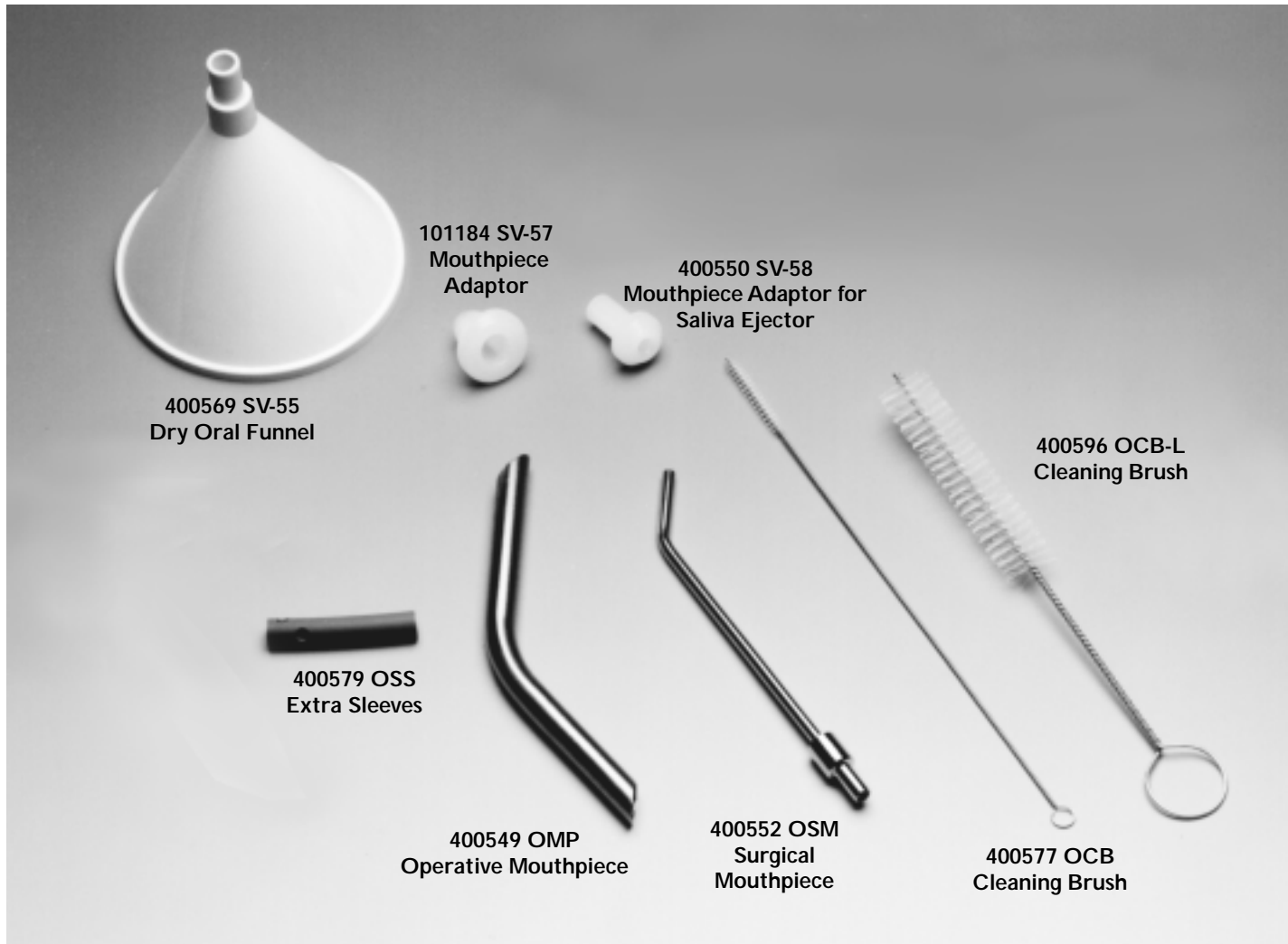

400569 SV-55 — Dry Oral Funnel - made of sterile white polycarbonate

400549 OMP — Operative Mouthpiece, stainless 7/16" O.D.

400552 OSM — Surgical Mouthpiece - fits standard and saliva ejector adaptors

400579 OSS — Extra Sleeves for Saliva Ejector (Packed 12)

400577 OCB — Mouthpiece Cleaning Brushes - for cleaning oral mouthpieces

400596 OCB-L — Mouthpiece Cleaning Brush - for cleaning high volume mouth pieces and valves. Brush diameter 3/4", brush length 8"

101184 SV-57 — Mouthpiece Adaptor, Standard

400550 SV-58 — Mouthpiece Adaptor for Saliva Ejector (Fits into SV-57 adaptor)

(800) 346-DENT

All prices are FOB Factory, Leicester, NY. Prices and designs are subject to change without notice.



FROM
ALLIED
DIAGNOSTIC IMAGING RESOURCES, INC.
MAKERS OF
REDI-CHEM

Steri-Dent • 2364 Leicester Road, P.O. Box 175 • Leicester, NY 14481-0175

Toll Free: (800) 346-DENT in US • Phone: (585) 382-3223 • Fax: (585) 382-3736 • email: kwhite@cpac-fuller.com